

Instruction Manual

CMS-50H



Instructions to User

Dear users, thank you very much for purchasing the Pulse Oximeter.

This Manual is written and compiled in accordance with the council directive MDD93/42/EEC for medical devices and harmonized standards. In case of modifications and software upgrades, the information contained in this document is subject to change without notice.

The Manual describes, in accordance with the Pulse Oximeter's features and requirements, main structure, functions, specifications, correct methods for transportation, installation, usage, operation, repair, maintenance and storage, etc. As well as the safety procedures to protect both the user and equipment. Refer to the respective chapters for details.

Please read the User Manual carefully before using this product. The User Manual which describes the operating procedures should be followed strictly. Failure to follow the User Manual may cause measuring abnormality, equipment damage and human injury. The manufacturer is NOT responsible for the safety, reliability and performance issues and any monitoring abnormality, human injury and equipment damage due to users' negligence of the operation instructions. The manufacturer's warranty service does not cover such faults.

Owing to the forthcoming renovation, the specific products you received may not be totally in accordance with the description of this User Manual. We would sincerely regret for that.

This product is medical device, which can be used repeatedly.

WARNING:

- **Discomfort or pain may appear if using the device ceaselessly, especially for the microcirculation barrier patients. It is recommended that the pulse oximeter should not be used on the same finger for more than 2 hours.**
- **For the special patients, there should be a more prudent inspecting in the placing process. The device can not be clipped on the edema and tender tissue.**
- **The light (the infrared is invisible) emitted from the device is harmful to the eyes, so the user and the maintenance man should not look at the light source.**
- **Testee can not use enamel or other makeup.**
- **Testee's fingernail can not be too long.**
- **Please refer to the correlative literature about the clinical restrictions and caution.**
- **This device is not intended for treatment.**

The User Manual is published by our company. All rights reserved.

Contents

| | | |
|-----------|---|-----------|
| 1 | Safety | 1 |
| 1.1 | Instructions for safe operations..... | 1 |
| 1.2 | Warning..... | 1 |
| 1.3 | Attention..... | 1 |
| 2 | Overview | 2 |
| 2.1 | Features..... | 2 |
| 2.2 | Major applications and scope of application..... | 2 |
| 2.3 | Environment requirements..... | 3 |
| 3 | Principle | 3 |
| 4 | Technical specifications | 3 |
| 4.1 | Main performance..... | 3 |
| 4.2 | Main Parameters..... | 4 |
| 5 | Installation | 4 |
| 5.1 | View of the front panel..... | 4 |
| 5.2 | Lanyard installation..... | 5 |
| 5.3 | Probe installation..... | 5 |
| 5.4 | Accessories..... | 5 |
| 6 | Operating Guide | 5 |
| 6.1 | Application method..... | 5 |
| 6.2 | Attention for operation..... | 10 |
| 6.3 | Clinical restrictions..... | 10 |
| 7 | Maintain, transportation and storage | 11 |
| 7.1 | Cleaning and Disinfecting..... | 11 |
| 7.2 | Maintain..... | 11 |
| 7.3 | Transportation and storage..... | 11 |
| 8 | Troubleshooting | 11 |
| 9 | Key of Symbols | 12 |
| 10 | Function Specification | 13 |
| | Appendix 1 | 14 |

1 Safety

1.1 Instructions for safe operations

- ✧ Check the main unit and all accessories periodically to make sure that there is no visible damage that may affect patient's safety and monitoring performance. It is recommended that the device should be inspected weekly at least. When there is obvious damage, stop using the device.
- ✧ Necessary maintenance must be performed by qualified service engineers ONLY. Users are not permitted to maintain it by themselves.
- ✧ The oximeter cannot be used together with devices not specified in User's Manual. Only the accessory that is appointed or recommendatory by manufacture can be used with this device.
- ✧ This product is calibrated before leaving factory.


1.2 Warning


- ⚠ Explosive hazard—DO NOT use the oximeter in the environment with tinder such as anesthetic .
- ⚠ DO NOT use the oximeter while the patient is being scanned by MRI or CT.
- ⚠ DO NOT strand the lanyard in order to avoid device drop and damage.The lanyard is made of non-sensitive material.Please do not use lanyard if the user is allergic to lanyard.Do not enwind neck with lanyard in order to avoid accident.
- ⚠ The person who is allergic to rubber can not use this device.
- ⚠ The disposal of scrap instrument and its accessories and packing (including battery, plastic bags, foams and paper boxes) should follow the local laws and regulations.
- ⚠ Please check the packing before use to make sure the device and accessories are totally in accordance with the packing list, or else the device may have the possibility of working abnormally.
- ⚠ Please choose the accessories which are appointed or recommended by the manufacturer for avoiding device damage.
- ⚠ Please don't measure this device with functional tester for the device's related information.


1.3 Attention


- 🔔 Keep the oximeter away from dust, vibration, corrosive substances, tinder, high temperature and moisture.
- 🔔 If the oximeter gets wet, please stop operation.
- 🔔 When it is carried from cold environment to warm or humid environment, please do not use it immediately.
- 🔔 DO NOT operate button on front panel with sharp things.
- 🔔 High temperature or high pressure steam disinfection for the oximeter is not permitted. Refer to User Manual in the relative chapter (7.1) for cleaning and disinfection.
- 🔔 Do not have the oximeter immerged in liquid. When it needs cleaning, please wipe its surface with medical alcohol by soft material. Do not spray any liquid on the device directly.
- 🔔 When cleaning the device with water, the temperature should be lower than 60°C.
- 🔔 The fingers which are too thin or too cold may affect the measure accuracy , please clip the thicker finger such as thumb and middle finger deeply enough into the probe.
- 🔔 The oximeter probe which is optional can be used to children.
- 🔔 The update period of data is less than 5 seconds, which is changeable according to different


individual pulse rate.


 Please read the measure value when the waveform on screen is equably and steady-going. This measure value is optimal value, and the waveform at the moment is the standard one.

 If some abnormal conditions appear on the screen during test process, pull the finger out and reinsert to restore normal use.

 The device has normal life for three years since the first electrified use.

 The device has alarm function, this function can either be paused, or closed (default setting) for good. Please check the chapter 6.1 as reference.

 The device has the function of beyond limit alarm. When the measure data is beyond the highest or lowest limit, the device would start alarm automatically on the premise of the alarming function is on.

 The device may be not fit for all patients. If you are unable to receive approving measure, discontinue use.

 Do not contort or drag the wire of the device.

2 Overview

The pulse oxygen saturation is the percentage of HbO₂ in the total Hb in the blood, so-called the O₂ concentration in the blood. It is an important bio-parameter for the respiration. A number of diseases relating to respiratory system may cause the decrease of SpO₂ in the blood, furthermore, some other causes such as the malfunction of human body's self-adjustment, damages during surgery, and the injuries caused by some medical checkup would also lead to the difficulty of oxygen supply in human body, and the corresponding symptoms would appear as a consequence, such as vertigo, impotence, vomit etc. Serious symptoms might bring danger to human's life. Therefore, prompt information of patients' SpO₂ is of great help for the doctor to discover the potential danger, and is of great importance in the clinical medical field.

The Pulse Oximeter is small in volume, low in power consumption, convenient in operation and portable. With high-definition display screen, the device is concise and fashion. It is only necessary for patient to put one finger into probe for diagnosis, and the display screen will directly show the SpO₂ value, pulse rate value, Perfusion Index value and pulse waveform with the high veracity and repetition.

2.1 Features

A Ultra-thin design, concise and fashion.

B Small in volume, light in weight and convenient in carrying.


C Low power consumption

D Display direction can be changed automatically, easy to view.

E Perfusion Index measure.

2.2 Major applications and scope of application

The Pulse Oximeter can be used in measuring the pulse oxygen saturation, pulse rate and perfusion index through finger. The product is fit for family, hospital, oxygen bar, community healthcare, physical care in sports (It can be used before or after doing sports and it is not recommended to use the device during the process of having sport) and etc.

 **The problem of overrating would emerge when the patient is suffering from toxicosis**

which caused by carbon monoxide, the device is not recommended to be used under this circumstance.

2.3 Environment requirements

Storage Environment

- a) Temperature : $-40^{\circ}\text{C}\sim+60^{\circ}\text{C}$
- b) Relative humidity : $5\%\sim95\%$
- c) Atmospheric pressure : $500\text{hPa}\sim1060\text{hPa}$

Operating Environment

- a) Temperature: $10^{\circ}\text{C}\sim40^{\circ}\text{C}$
- b) Relative Humidity : $30\%\sim75\%$
- c) Atmospheric pressure: $700\text{hPa}\sim1060\text{hPa}$

3 Principle

Principle of the Oximeter is as follows: An experience formula of data process is established taking use of Lambert Beer Law according to Spectrum Absorption Characteristics of Reductive Hemoglobin (Hb) and Oxyhemoglobin (HbO_2) in glow & near-infrared zones. Operation principle of the device is: Photoelectric Oxyhemoglobin Inspection Technology is adopted in accordance with Capacity Pulse Scanning & Recording Technology, so that two beams of different wavelength of lights can be focused onto human nail tip through perspective clamp finger-type sensor. Then measured signal can be obtained by a photosensitive element, information acquired through which will be shown on screen through treatment in electronic circuits and microprocessor.

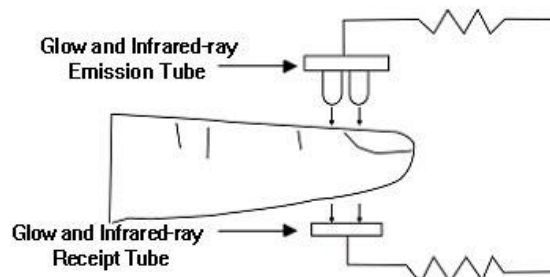


Figure 1 work elements

4 Technical specifications

4.1 Main performance

- SpO_2 value display
- Pulse rate value display, bar graph display
- Pulse waveform display
- Perfusion Index value display
- Charge function: Lithium battery is the power supply and can be charged up time and again.
- Low-power indication: low-power indication symbol appears before working abnormally which is due to low-power.
- Automatic power-off when there is no finger in device under measure interface in 5 seconds.
- With direction sensor, the display direction can be changed by automatic or manual mode
- Screen brightness can be adjusted
- Pulse sound indication

- With alarm function
- With SpO₂ value and pulse rate value record function
- The record data can be uploaded to the computer
- It can be connected with oximeter probe

4.2 Main Parameters

A Measurement of SpO₂

Measuring range: 0% ~ 100%

Accuracy: 70~100%:±2%;Below 70%:unspecified.

B Measurement of pulse rate

Measuring range:30bpm~250bpm

Accuracy: 30bpm~100bpm:±2bpm;101bpm~250bpm: ±2%

C Perfusion Index

Range: 0% ~ 20%

D Resolution

SpO₂ : 1%, Pulse rate: 1bpm

E Measurement Performance in Weak Filling Condition:

SpO₂ and pulse rate can be shown correctly when pulse-filling ratio is 0.4%. SpO₂ error is ±4%; pulse rate error is ±2 bpm or ±2% (select larger).

F Resistance to surrounding light:

The deviation between the value measured in the condition of man-made light or indoor natural light and that of darkroom is less than ±1%.

G Power supply requirement: 3.6 V DC ~4.2V DC.

H Optical Sensor

Red light (wavelength is 660nm, 6.65mW)

Infrared (wavelength is 880nm, 6.75mW)

I Adjustable alarm range:

SpO₂ : 0%~100%

Pulse Rate: 0bpm~254bpm

5 Installation

5.1 View of the front panel

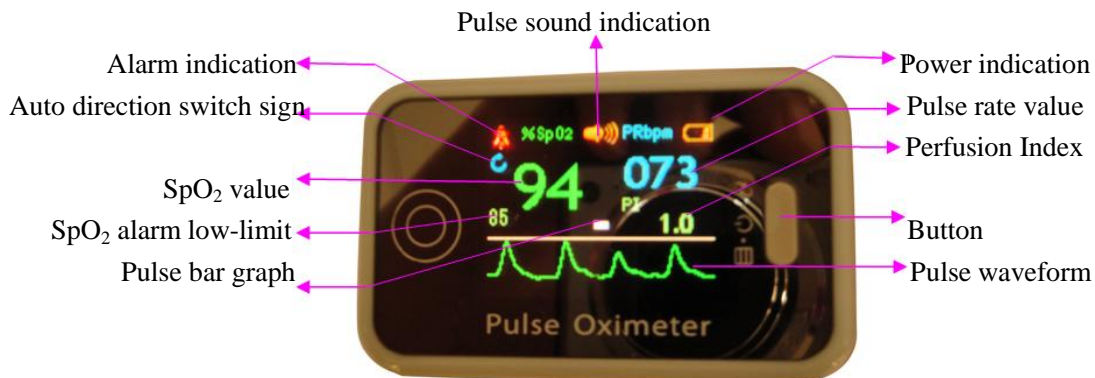



Figure 2. Front View

5.2 Lanyard installation

- A. Put the thinner side of the lanyard through the hole.
- B. Put the wider side of the lanyard through the thinner side which has been put through the hole, then tighten it.

5.3 Probe installation

- A. Open the USB cover,inserting the probe into the USB port of the pulse oximeter.
- B. Insert the finger into the probe (the side with finger sign and the nail should be in the same side).The connected probe could work normally when the luminescent tube in the device doesn't emit red light.

 In order to ensure device life and measure accuracy,external probe is limited to the matching model ,and can't be replaced by other model.Or else the device will be damaged or usage is affected.

5.4 Accessories

- A. A lanyard
- B. A User Manual
- C. A power adapter(optional)
- D. A data line
- E. A disk (PC software)
- F. A oximeter probe (optional)

6 Operating Guide

6.1 Application method

6.1.1 Measurement

- a. Insert the finger into the probe of the device.
- b. Long press button to turn the device on.
- c. Do not shake the finger and keep the patient in a stable state during the process.
- d. The data can be read directly from the display screen in the measuring interface.

6.1.2 Lay finger

The right method of laying finger is as figure 3.



Figure 3

⚠️Fingernails and the luminescent tube should be in the same side.

⚠️Automatic power-off when there is no finger in device in 5 seconds.

6.1.3 Pause alarm

- a. Alarm includes the alarm of measure data's going beyond the limits, the alarm of low-power, and the alarm of finger out.
- b. In the measuring interface, if the alarm function is on, during the period of alarming, alarm can be suspended by short pressing the button, but the function will be renewed in about 60 seconds.
- c. If you want to turn off the alarm, you should enter the menu for operation. Please refer to chapter "Alarm setting" for detail.

6.1.4 Change display direction

The device could change display direction by automatic or manual mode. When the direction sensor is on, the device could change display direction according to the handing direction; When the direction sensor is off, short press button to change the display direction. Please refer to chapter "Direction sensor setting" for detail.

6.1.5 Menu operations

In the measuring interface, long press button to enter the main menu interface as figure 4. (When the display direction is lengthways, you can not enter the main menu interface. Please change to landscape orientation.). The user can adjust the settings through the main menu, such as backlight, alarm, direction sensor, data transmission (with the data line). The specific operation methods are as follows:

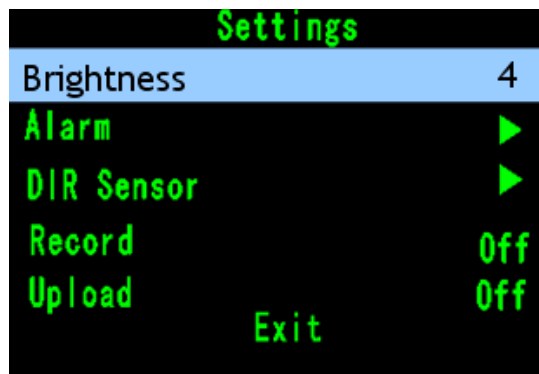


Figure 4 Main Menu Interface

A Brightness adjustment

In the main menu interface, short press button to move the choice bar to the "Brightness" item, then long press the button to adjust the brightness of screen (4 levels adjustable).

B Alarm setting

In the main menu interface, short press button to move the choice bar to "Alarm" item, then long press button to enter the alarm setting menu as figure 5:

| Dir | Up |
|-------------------------|-----|
| SpO ₂ ALM HI | 100 |
| SpO ₂ ALM LO | 85 |
| PR ALM HI | 120 |
| PR ALM LO | 50 |
| Alarm | Off |
| Pulse Sound | Off |
| Exit | |

Figure 5. Alarm Setting Menu

a The high/low limit of alarm setting

In alarm setting menu, short press button to move the choice bar to "Dir" item, long press button to choose Up or Down (This will be the direction the value of the high-low limits of SpO₂ and pulse rate will be adjusted). Short press button to adjust alarm item: SpO₂ high limit (SpO₂ ALM HI), SpO₂ low limit (SpO₂ ALM LO), Pulse rate high limit (PR ALM HI), Pulse rate low limit (PR ALM LO). Long press button to change the value. Long press button once, and the value raise or descend once. The low limit can't be above high limit, and the high limit can't be below the low limit.

⚠️ If the alarm function is on, the device will provide medium-priority alarm signal when the data of SpO₂ or pulse rate is beyond the limit. Intermittent alarm will occur and the values show in yellow.

Medium priority indicating that prompt operator response is required.

b The alarm indication setting

In alarm setting menu, short press button to move the choice bar to the "Alarm" item, then long press button to turn on or turn off the alarm sound. Choose "on" to turn on the alarm, and choose "off" to turn off the alarm.

c Pulse sound indication setting

In alarm setting menu, short press button to move the choice bar to the "Pulse Sound" item, then long press button to turn on or turn off the pulse sound. Choose "on" to turn on pulse sound, and choose "off" to turn off pulse sound.

d Exit operation

In alarm setting menu, short press button to move the choice bar to "Exit" item, then long press button to exit alarm setting menu and return to the main menu.

C Direction sensor setting

In the main menu interface, short press button to move the choice bar to "DIR Sensor" item, then long press button to enter the direction sensor setting menu as figure 6.

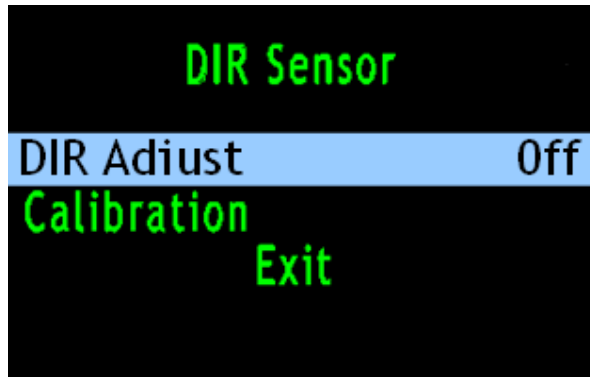


Figure 6 Direction sensor setting menu

a Turn on/off the direction sensor

In the direction sensor setting interface, short press button to move the choice bar to "DIR Adjust" item. Long press button to adjust this function. When you choose "auto" item, there is a auto direction switch sign in the measuring interface and the device could adjust display direction according to the handing direction. When you choose "off" item, in the measuring interface, auto direction switch sign disappears, and short press button to change the display direction.

b Direction sensor calibration function

When the direction sensor error is larger, the user could use this function to calibrate direction sensor.

In the sensor setting interface, short press button to move the choice bar to "calibration" item. Long press button, the display screen will appear prompt "You will calibrate the direction sensor". If you want to calibrate, please put the device horizontally, then move the choice bar to "yes" item and long press button to affirm. Wait 5 seconds until the display screen appears "Calibrate successfully", here calibrate direction sensor successfully. If prompt "Calibrate unsuccessfully" appears, this calibration is unsuccessful. The user could calibrate again according to the above operation until calibration is successful.



When calibrating, please put the device horizontally and stop measuring.

D Data record setting

This device can record 24 hours data including pulse rate and SpO₂ value accurately when the battery is full and upload the data to the computer with data line for display and analysis.

a In the main menu interface, short press button to move the choice bar to "Record" item, then long press button to enter the record time setting dialog box as figure 7.

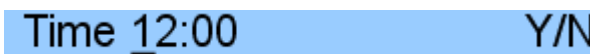


Figure 7 Record time setting dialog box

b Short press button to move the underline to the number that you want to set, then long press button to set time. After setting time, move the underline to "Y", then long press button to set "record" item as "on" and begin to record. If move the underline to "N", then long press button to cancel record. When recording, long press button to set "record" item as "off", here record function is closed, data is saved.

- c If the record function is turned on, a flashing red dot would appear on the screen when returning to the measure interface ,which means the device is recording.
- d If the device is recording,whatever interface the device is in (measuring interface, menu interface), the sign "Recording" would appear on the screen in 30 seconds, then the screen will be automatically shut down. If short press button at this moment,the sign "Recording" would appear on the screen,then the screen will be automatically shut down again; if long press the button, the device would return to the former interface.
- e If turning on the record function, the former saved data will be automatically deleted.
- f When recording, the pulse sound indication would be turned off for saving power, after the screen is shut down automatically.
- g. When the storage space is full, it displays “Memory is full” on the screen, and then shut down in a few seconds. But it will display “Memory is full” by the next time you turn on the device on the purpose of warning the user.

E Data transmission setting

Please install the affiliated software into the computer, and two icons would appear on the desktop after installation. The icon of "SpO₂ "is a program for receiving real-time data which is shown as figure 8; the icon of "SpO₂ Review" is a program for receiving stored data which is shown as figure 9.



Figure 8 SpO₂ program

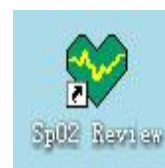


Figure 9 SpO₂ Review program

a Real-time data transmission

1. Please connect the device to computer with the data line which is affiliated with the device, then double click "SpO₂ "icon to open"SpO₂ " program.
2. On the display of the computer,you could see the data information of SpO₂ and pulse rate.
3. When the data line is pulled out from the computer,there will be a "Save data at view" dialog box on the display of the computer,and you can input the patient information in it.

b Record data transmission

- 1.Please connect the device to computer with the data line which is affiliated with the device, then double click "SpO₂ Review "icon to open"SpO₂ Review" program.Click the ‘New Session’ Icon in the software, and input the patient data ,then click ‘ok’.Here the Software will display “device connected, waiting for data”.
2. Long press button to enter the “Main Menu”and short press button to select “Upload” item, then long press button to upload the record data to computer for display and analysis.
3. When recording,, it is not applicable for the user to upload the storage data to computer.
4. After accomplishing the uploading,the information of SpO₂ and pulse rate recorded can be displayed on the display of the computer.

⚠️If the user choose to turn on the computer display function , it would probably take several

seconds for the data to appear on the computer screen.


F Exit the main menu


In the main menu interface, move the choice bar to "Exit" item, then long press button to exit the main menu and return to the measure interface.

6.1.6 Charge

There are two kinds of charge method:

- a) Connect the device to computer with data line, then the device should be in charge state.
- b) Connect the device to power supply with power adaptor, then the device should be in charge state.

 **In the charge state, the blue indication light shining means the device is charging up, the blue indication light quenching means the charge has been accomplished.**

 **If the alarm function is on, the device will provide high-priority alarm signal when the battery is in low power status. Intermittent alarm will occur. High priority indicating requires that operator responds immediately.**

6.2 Attention for operation

- A. Please check the device before using, and confirm that it can work normally.
- B. The finger should be in a proper position (see the attached illustration as figure 3 for reference), or else it may result in inaccurate measure.
- C. The ray between luminescent tube and photoelectric receiving tube must get across subject's arteriole.
- D. The oximeter should not be used at a location or limb tied with arterial canal or blood pressure cuff or receiving intravenous injection.
- E. Ensure nothing, such as a plaster, can impede the light passage., or else it may result in inaccurate measure of SpO₂, pulse rate and PI..
- F. Excessive ambient light may affect measurement accuracy. It includes fluorescent lamp, dual ruby light, infrared heater, direct sunlight and etc.
- G. Intense activity of the subject or extreme electrosurgical interference may also affect the accuracy.
- H. Testee can not use enamel or other makeup.
- I. Please clean and disinfect the device after operating according to the User Manual (7.1).

6.3 Clinical restrictions

- A. As the measure is taken on the basis of arteriole pulse, substantial pulsating blood flow of subject is required. For a subject with weak pulse due to shock, low ambient/body temperature, major bleeding, or use of vascular contracting drug, the SpO₂ waveform (PLETH) will decrease. In this case, the measurement will be more sensitive to interference.
- B. For those with a substantial amount of staining dilution drug (such as methylene blue, indigo green and acid indigo blue), or carbon monoxide hemoglobin (COHb), or methionine (Me+Hb) or thiosalicylic hemoglobin, and some with icterus problem, the SpO₂ determination by this device may

be inaccurate.

C. The drugs like dopamine, procaine, prilocaine, lidocaine and butacaine may also be a major factor resulted in serious error of SpO₂ measure.

D. The SpO₂ value serves as a reference value for judgment of anemic anoxia and toxic anoxia, for some patients with serious anemia may also report good SpO₂ measurement.


7 Maintain、 transportation and storage

7.1 Cleaning and Disinfecting

Using medical alcohol to wipe the device for disinfecting, nature dry or clean it with clean soft cloth.

7.2 Maintain

A Please clean and disinfect the device before using according to the User Manual (7.1).

B Please recharge the battery when the screen shows low-power (the battery power is ).

C Recharge the battery soon after the over-discharge. The device should be recharged every six months when it is not regular used. It can extend the battery life following this guidance.

D The device needs to be calibrated once a year (or according to the calibrating program of hospital). It also can be performed at the state-appointed agent or just contact us for calibration.

7.3 Transportation and storage

A. The packed device can be transported by ordinary conveyance or according to transport contract. The device can not be transported mixed with toxic, harmful, corrosive material.











B. The packed device should be stored in room with no corrosive gases and good ventilation. Temperature: -40°C~60°C; Relative Humidity: ≤95%


8 Troubleshooting

| Trouble | Possible Reason | Solution |
|--|--|---|
| The SpO₂ and Pulse Rate can not be displayed normally. | The finger is not properly positioned. | Place the finger properly and try again. |
| | The patient's SpO ₂ is too low to be detected. | Try again; Go to a hospital for a diagnosis if you are sure the device works all right. |
| The SpO₂ and Pulse Rate are not displayed stably. | The finger is not placed inside deep enough. | Place the finger properly and try again. |
| | The finger is shaking or the patient is moving. | Let the patient keep calm. |
| The device can not be turned on. | The battery is drained away or almost drained away. | Please charge up battery. |
| | The malfunction of the device. | Please contact the local service center. |
| The display is off suddenly | The device is set to shut down automatically in 5 seconds when there is no signal. | Normal. |

| | | |
|---|---|--|
| | The battery is drained away or almost drained away. | Please charge up battery. |
| The device can not be used for full time after charge | The battery is not full charged. | Please recharge the battery |
| | The battery is broken | Please contact the local service center. |
| The battery can not be full charged even after 10 hours charging time. | The battery is broken | Please contact the local service center. |

9 Key of Symbols

| Signal | Description |
|---|--|
|  | Warning – See User Manual |
| %SpO2 | The pulse oxygen saturation (%) |
| PRbpm | Pulse rate (bpm) |
| PI | Perfusion Index (%) |
|  | The battery power is full |
|  | Low-power |
|  | Close the alarm sound indication |
|  | Pause the alarm sound indication |
|  | Open the alarm sound indication |
|  | Close the pulse sound indication |
|  | Open the pulse sound indication |
|  | Auto direction switch sign |
|  | Menu button/power button/function button |

| | |
|---|-------------------------|
| IPX1 | Ingress of liquids rank |
|  | USB |

10 Function Specification

| Information | Display Mode |
|--|--|
| The Pulse Oxygen Saturation (SpO ₂) | 2-digit digital OLED display |
| Pulse Rate (PR) | 3-digit digital OLED display |
| Perfusion Index (PI) | 3-digit digital OLED display |
| Pulse Intensity (bar-graph) | bar-graph OLED display |
| SpO₂ Parameter Specification | |
| Measuring range | 0% ~ 100%, (the resolution is 1%). |
| Accuracy | 70% ~ 100%: ±2% ,Below 70% unspecified. |
| Pulse Parameter Specification | |
| Measuring range | 30bpm ~ 250bpm, (the resolution is 1bpm) |
| Accuracy | 30bpm ~ 100bpm: ±2bpm; 101bpm ~ 250bpm: ±2% |
| Perfusion Index Specification | |
| Range | 0% ~ 20%. |
| Safety Type | Interior Battery, BF Type |
| Pulse Intensity | |
| Range | Continuous bar-graph display, the higher display indicates the stronger pulse. |
| Battery Requirement | |
| Voltage 3.7 rechargeable lithium battery × 1 (The red wire on the battery denotes anode, the black wire on the battery denotes cathode.) | |

| | |
|--|------------------------------------|
| Battery working life | |
| Charge and discharge no less than 500 times. | |
| Power Adapter (selected) | |
| Output voltage | DC 5V |
| Output current | 250mA |
| Dimensions and Weight | |
| Dimensions | 58(L) × 36(W) × 26(H) mm |
| Weight | About 45g (with a lithium battery) |

Appendix 1

| State | Alarm condition delay | Alarm signal generation delay |
|------------------------|-----------------------|-------------------------------|
| Low-power alarm | 1s | 20ms |
| Spo ₂ alarm | 330ms | 20ms |
| Pulse rate alarm | 330ms | 20ms |
| Probe error alarm | 16ms | 20ms |